

# Reusable EEG Cup Electrode Silver/Silver Chloride



**Technomed understands the importance of a crystal clear signal, ease of use and electrode durability.**

**Technomed's Reusable Ag/AgCl EEG Cup Electrode is designed to provide a reliable signal and longtime ease of use.**

**The electrode head is designed with a larger dome volume for long term use.**

**The connector is optimized for easy headbox insertion and extraction.**


## **User Benefits**

- Clear and reliable signal with optimized layer of frequency stable Ag/AgCl on a solid Silver cup.
- Extra large dome volume for longer use compared to a standard dome.
- Short cup-to-cable connection to reduce lift risk.
- Optimal patient comfort with soft and flexible lead wires.
- Optimized connector for easy headbox insertion and extraction.
- Long electrode life with watertight design and Kevlar reinforced wires.
- Available in two different cup sizes.

# Available variants

Materials	Cup Size	Cup Height	Wire Length	Packing	Part Number
Silver/Silver Chloride (Ag/AgCl)	10 mm	3.3 mm	1 m	Sets of 12, 2 x 6 colors	TE/C22-434/99
Silver/Silver Chloride (Ag/AgCl)	10 mm	3.3 mm	1.5 m	Sets of 12, 2 x 6 colors	TE/C22-634/99
Silver/Silver Chloride (Ag/AgCl)	10 mm	3.3 mm	2 m	Sets of 12, 2 x 6 colors	TE/C22-834/99
Silver/Silver Chloride (Ag/AgCl)	10 mm	3.3 mm	2.5 m	Sets of 12, 2 x 6 colors	TE/C22-934/99
Silver/Silver Chloride (Ag/AgCl)	6 mm	2.3 mm	1.5 m	Sets of 12, 2 x 6 colors	TE/C24-634/99

**2 x 6 colors:** red, green, white, yellow, blue, brown

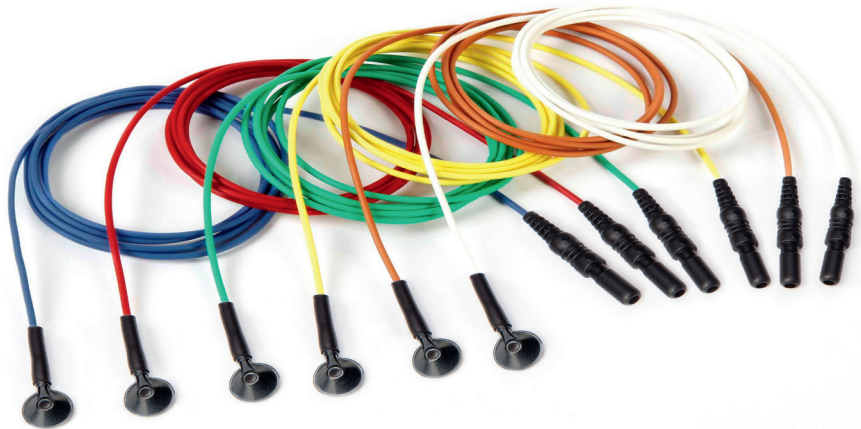

## Specifications

Cup	Sterling Silver / AgCl layer
Not made with natural rubber latex	√
REACH & RoHS (EU) compliant	√
Lead wire design	Tinsel wire
Lead wire core	Kevlar and Silver Plated Copper
Lead wire insulation	TPU
Spring loaded connector	Gold Plated Brass / PVC
Connector	DIN 42802-1

**Cleaning and disinfection:** see IFU for details

Your Technomed Business Partner:



[www.technomed.nl](http://www.technomed.nl)

Legal manufacturer:

Technomed Europe

Amerikalaan 71

6199 AE | Maastricht-Airport | The Netherlands

T +31 43 408 68 68 | [info@technomed.nl](mailto:info@technomed.nl)

Technomed USA

17200 Medina Road, Suite 800

Plymouth | MN 55447 | United States

T +1 651 288 0890 | [info@technomed-usa.com](mailto:info@technomed-usa.com)